


DimAir Line

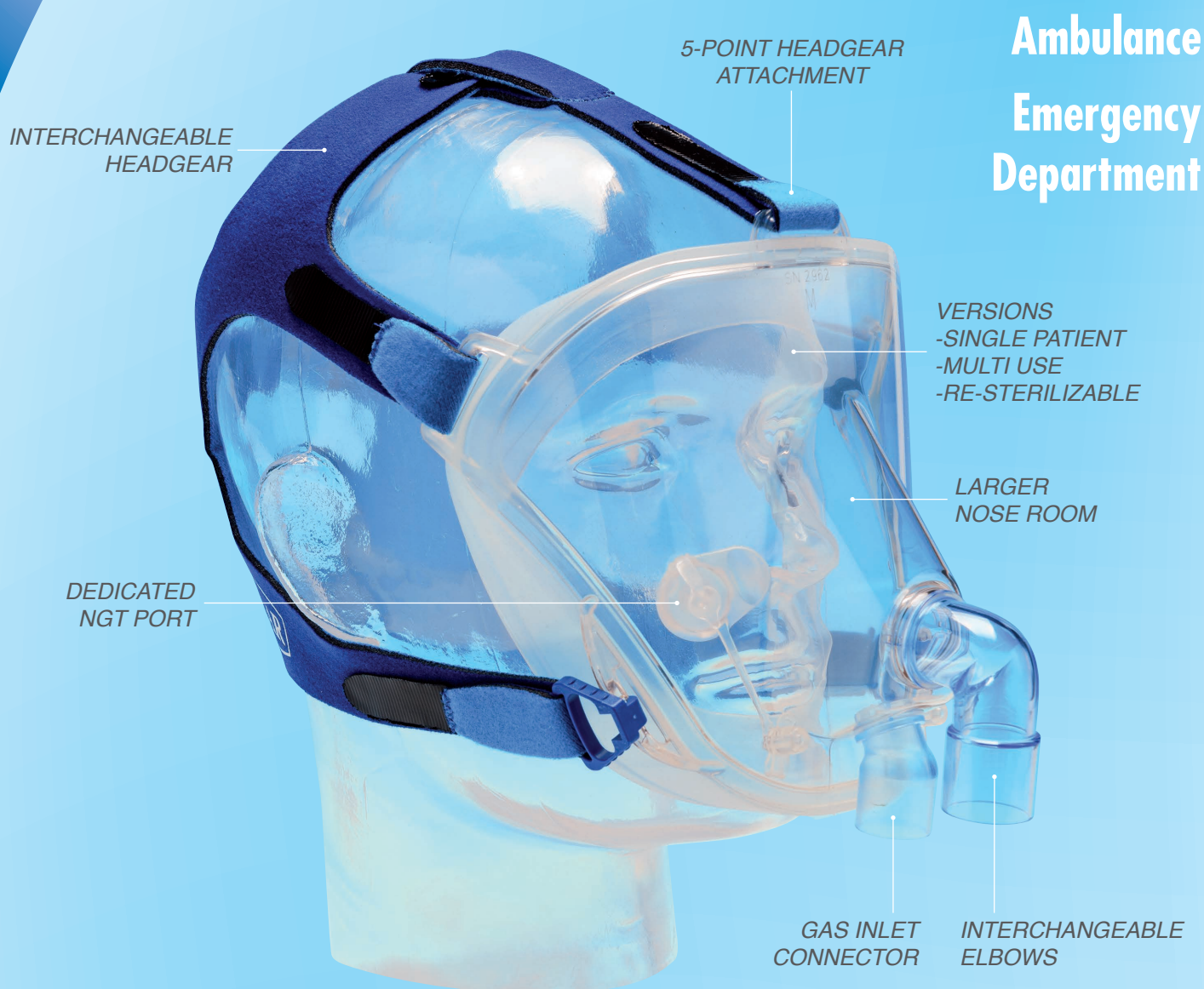


TOTAL FACE
DIMAX ZERO

A new generation total face mask for Non-Invasive Ventilation in PSV, BIPAP and CPAP mode. Its original configuration reduces to zero the CO₂ rebreathing. Available in different sizes to fit any paediatric or adult patient, DiMax can also be used in combination with dual or single limb ventilators as well as CPAP high flow generators by simply changing the elbow. Unique features such as the 5-point mask attachment, NGT airtight port, bronchoscopy airtight port, extra nose room, make DiMax Zero the perfect interface for Non Invasive Ventilation.

8 plus

- zero rebreathing
- securement points
- comfort
- sealing
- practicality
- versatility
- safety
- versions

Fluid Dynamics of CO₂ rebreathing in standard total face mask compared to DiMax Zero

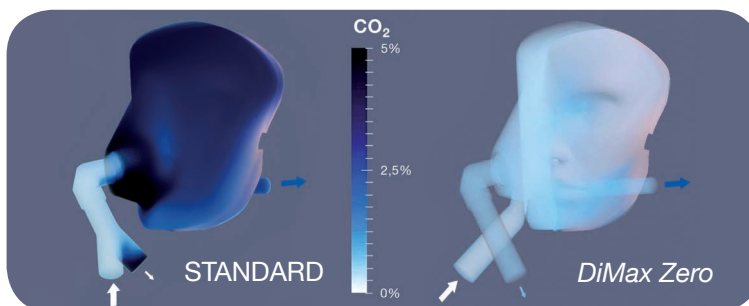

TECHNICAL FEATURES

- Zero dead space
- 5-point headgear attachment
- Larger nose room
- Dedicated ports for NG tube and bronchoscopy probe
- Single patient and reprocessable versions
- Interchangeable elbows:
 - 22F non-vented (blue)
 - 22M vented (clear)

AVANTAGES

- Zero CO₂ rebreathing
- Improved patient/ventilator synchronism
- Improved airtight performance
- Improved comfort regardless of patient face contour or nose size
- Less headgear pressure
- Less risk of skin lesions
- NG tubes and bronchoscopy probes can be used without air leaks or disrupting the ventilation
- Single or multi patient versions, reprocessable in autoclave @ 134° - Sterrad - ETO in autoclave @ 134° - gas plasma - ETO



The fifth point of attachment allows a better sealing onto the patient's forehead



Dedicated air tight NG tube port



"Non-Vented" elbow for single limb ventilator



Configuration for CPAP with "Easy Vent" Venturi Generator



Configuration with dual limb ventilator



Configuration with single limb ventilator or CPAP Generator

REFERENCES

CODES*			SIZE
SINGLE PATIENT**	MULTI PATIENT**	REPROCESSABLE***	
700/6797	700/6797P	700/6801	SMALL
700/6798	700/6798P	700/6802	MEDIUM
700/6799	700/6799P	700/6803	LARGE
700/6800	700/6800P	700/6804	EXTRA LARGE

* Including "non-vented" blue elbow, "vented" clear elbow and 1 headgear
 ** single packed 5 pcs multi box
 *** single packed